SUPER ERECTA SHELF®

Inhalation Therapy Cart

☐ Super Erecta Shelf® Construction:
Allows custom height adjustment and positive locking without the use of nuts, bolts or setscrews and without the use of special tools.


Note: To maintain the conductivity of the cart, aluminum split sleeves must be used rather than plastic.

- □ PVC-Coated Top & Bottom Inserts: The compartmented inserts, which snap over the top and bottom frames, have a resilient PVC coating which aids in reducing vibration and noise.
- ☐ Chrome-Plated Finish on handles, frames and shelf.
- □ Accommodates Any Combination of twelve D & E cylinders. Each PVC-coated insert has twelve compartments, each measuring 4¹/2" (114mm) square.
- 5" (127mm) Conductive Casters. See caster NOTICE on assembly instruction sheet.

Copyright © 1984 InterMetro Industries Corp.



ITC12C

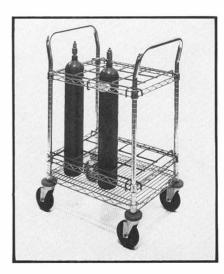



30.40

DIMENSIONS

	Width		Length		Height		Approx. Pkd. Wt.	
Cat. No.	(In.)	(mm)	(in.)	(mm)	(in.)	(mm)	(lbs.)	(kg)
ITC12C	201/2	520	26	660	39	990	401/2	18

